

# ethidium monoazide bromide (EMA)

## Description

**Ethidium monoazide** is a fluorescent photoaffinity label that, after photolysis, binds covalently to nucleic acids both in solution and in cells that have compromised membranes. The fluorescence of ethidium monoazide is weak, but the intensity increases ~15-fold on binding to DNA with excitation/emission maxima of ~504/600 nm.

Catalog Number  
E1374

Size  
5 mg

List Price  
(USD) 189.00

## Specifications

### General Specifications

Color:	Orange	
Platform:	Flow Cytometer	
Product Size:	5 mg	
Detection Method:	Fluorescent	/
Cell Permeability:	Cell-Impermeant	/
Sub-Cellular Localization:	Nucleic Acids, Nucleus	



## Chemical Structures

ethidium monoazide bromide

Molecular Formula:


$C_{21}H_{18}BrN_5$

Molecular Weight:

420.3107

CAS Name/Number:

Phenanthridium, 3-amino-8-azido-5-ethyl-6-phenyl, bromide/ 58880-05-0

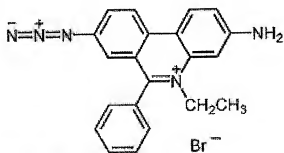

[Enlarge](#)

---

## Documents

Material Safety Data Sheets (MSDS)

- [E1374](#)

Certificates of Analysis (COA)

Molecular Probes Handbook

- [Photoactivatable Reagents, Including Photoreactive Crosslinkers and Caged Probes—Section 5](#)
- [Nucleic Acid Stains—Section 8.1](#)
- [Viability and Cytotoxicity Assay Reagents—Section 15.2](#)

## Additional Information



## Citations & References (96)

- Novel pentablock copolymers for selective gene delivery to cancer cells.

### Authors

Zhang B, Kanapathipillai M, Bisso P, Mallapragada S,

### Journal

Pharm Res (2009) 26:700-713

### Product Usage

Labeling of poly(diethylaminoethylmethacrylate)/Pluronic F127 copolymers for visualization of microscopy.

### ID:

PN67213

- Analysis of HCV-specific T cells by flow cytometry.

### Authors

Shiina M, Rehmann B,

### Journal

Methods Mol Biol (2009) 510:415-426

### Product Usage

Flow cytometric analysis of HCV-specific T-cell proliferation.

### ID:

PN67409

- In vitro micronucleus assay scored by flow cytometry provides a comprehensive evaluation

### Authors

Bryce SM, Bemis JC, Avlasevich SL, Dertinger SD

### Journal

Mutat Res (2007) 630:78-91

### ID:

PN58978

- Cationic albumin-conjugated pegylated nanoparticles allow gene delivery into brain tumors

### Authors

Lu W, Sun Q, Wan J, She Z, Jiang XG

### Journal

Cancer Res (2006) 66:11878-11887

### ID:

PN58479

- Selective removal of DNA from dead cells of mixed bacterial communities by use of ethidium bromide

### Authors

Nocker A, Camper AK

### Journal

Appl Environ Microbiol (2006) 72:1997-2004

### ID:

PN59007



## Product Reviews (0)

## Help Center



### Order Support

Search frequently asked questions, register as a customer, change your account settings, track your orders, report a problem, or simply find out how to contact us.



### Technical Support

Tools and resources for the Support team. Find product guides, troubleshooting

People Who Viewed

This Item Also Viewed

Image-iT® DEAD Green™ viability stain \*1 mM solution in DMSO\*

Cat #

I10291

Size

1 vial

List Price

**263.00**

(USD)



LIVE/DEAD® Reduced Biohazard Cell Viability Kit #1 \*green and red fluorescence\* \*100 assays\*

Cat #

L7013

Size

1 kit

List Price

**209.00**

(USD)



SYTOX® Green nucleic acid stain - 5 mM solution in DMSO

Cat #

S7020

Size

250 µl

List Price

**184.00**

(USD)

## Product Categories

- Nucleic Acid Cell Stains




## **Related Applications**

- [Cell Structure](#)
- [Cell Tracing & Tracking](#)
- [Cell Viability, Proliferation & Function](#)

## **Related Products**

- [Nuclear Probes](#)
- [Structural Probes](#)



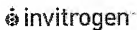

# nuclear yellow (Hoechst S769121, trihydrochloride, trihydrate)

## Description

**Nuclear yellow (Hoechst S769121)** exhibits excitation/emission maxima ~335/495 nm when bound to DNA.

Catalog Number  
N21485

Size  
10 mg

List Price  
(USD) 105.00

---

## Specifications

### General Specifications

Color:	Yellow
Platform:	Fluorescence Microscope
Product Size:	10 mg
Detection Method:	Fluorescent /
Cell Permeability:	Cell-Impermeant /
Sub-Cellular Localization:	Nucleic Acids, Nucleus



## Chemical Structures

nuclear yellow

Molecular Formula:


$C_{25}H_{34}Cl_3N_7O_5S$

Molecular Weight:

651.01

CAS Name/Number:

Benzenesulfonamide, 4-[5-(4-methyl-1-piperazinyl)[2,5'-bi-1H-benzimidazol]-2'-yl]-, trihydrochlorid

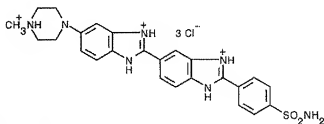

Q. Enlist

---

## Documents

### Manuals & Protocols

- [Hoechst Stains](#)

### Material Safety Data Sheets (MSDS)

- [nuclear yellow \(Hoechst S769121, trihydrochloride, trihydrate\)](#)

### Certificates of Analysis (COA)

### Molecular Probes Handbook

- [Probes for the Nucleus—Section 12.5](#)
- [Nucleic Acid Stains—Section 8.1](#)

## Additional Information



## Citations & References (189)

- Specific heterochromatic banding of metaphase chromosomes using nuclear yellow.

### Authors

Pinna-Senn E, Lisanti JA, Ortiz MI, Dalmasso G, Bella JL, Gosálvez J, Stockert JC

### Journal

Biotech Histochem (2000) 75:132-140

### ID:

PN38454

- The limbic zone of the rabbit and rat claustrum: a study of the claustrorostromal connective transport of fluorescent tracers.

### Authors

Majak K, Kowiński P, Morýs J, Spodnik J, Karwacki Z, Wisniewski HM

### Journal

Anat Embryol (Berl) (2000) 201:15-25

### ID:

PN39383

- Bifurcating projections from the cerebellar fastigial neurons to the thalamic supragenicular nuclei.

### Authors

Katoh YY, Arai R, Benedek G

### Journal

Brain Res (2000) 864:308-311

### ID:

PN41558

- Oncogenic Ras induces p19ARF and growth arrest in mouse embryo fibroblasts lacking p16INK4a.

### Authors

Groth A, Weber JD, Willumsen BM, Sherr CJ, Roussel MF

### Journal

J Biol Chem (2000) 275:27473-27480

### ID:

PN43117

- Higher activities of acetylcholinesterase and choline acetyltransferase in jaw-opening muscles in the rat.

### Authors

Kawagishi S

### Journal

Arch Oral Biol (1999) 44:197-200

### ID:

PN39382

## Product Reviews (0)



## Help Center



### Order Support

Search frequently asked questions, register as a customer, change your account settings, track your orders, report a problem, or simply find out how to contact us.



### Technical Support

Tools and resources for the Support team. Find product guides, troubleshooting

People Who Viewed

This Item Also Viewed

Hoechst 33342, trihydrochloride, trihydrate - 10 mg/mL solution in water

Cat #

H3570

Size

10 ml

List Price

79.00

(USD)



Image-iT® DEAD Green™ viability stain \*1 mM solution in DMSO\*

Cat #

I10291

Size

1 vial

List Price

263.00

(USD)



TO-PRO®-3 iodide (642/661) - 1 mM solution in DMSO

Cat #

T3605

Size

1 ml

List Price

255.00

(USD)

## Product Categories


- Nucleic Acid Cell Stains



## Related Applications

- [Cell Structure](#)
- [Cell Tracing & Tracking](#)



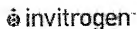

# propidium iodide

## Description

Propidium iodide is a popular red-fluorescent nuclear and chromosome counterstain. Since propidium iodide is not permeant to live cells, it is also commonly used to detect dead cells in a population. This dye is also available in solution (P-3566).

Catalog Number

P1304MP

Size

100 mg



List Price

(USD) 119.00

---

## Specifications

### General Specifications

Color:	Red
Platform:	Flow Cytometer, Fluorescence Microscope
Product Size:	100 mg
Detection Method:	Fluorescent 
Cell Permeability:	Cell-Impermeant 
Sub-Cellular Localization:	Cytoplasm & Cytosol



## Chemical Structures

propidium iodide

Molecular Formula:


$C_{27}H_{34}I_2N_4$

Molecular Weight:

668.4

CAS Name/Number:

Phenanthridinium, 3,8-diamino-5-[3-(diethylmethylammonio)propyl]-6-phenyl-, diiodide/ 25535-16-

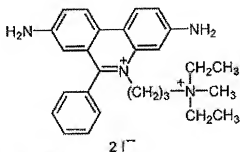


[Enlarge](#)

## Fluorescence Spectra

Propidium iodide/DNA

Absorption and fluorescence emission spectra of propidium iodide bound to DNA.

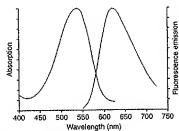

[Enlarge](#)



## Documents

### Manuals & Protocols

- [Propidium Iodide Nucleic Acid Stain](#)
- [Alexa Fluor® 488 Annexin V/Dead Cell Apoptosis Kit](#)

### Material Safety Data Sheets (MSDS)

- [PI304MP](#)

### Certificates of Analysis (COA)

### Molecular Probes Handbook

- [Probes for the Nucleus—Section 12.5](#)
- [Nucleic Acid Stains—Section 8.1](#)
- [Polar Tracers—Section 14.3](#)
- [Viability and Cytotoxicity Assay Reagents—Section 15.2](#)

## Additional Information

### Citations & References (1120)

- [Lack of p21 expression links cell cycle control and appendage regeneration in mice.](#)

#### Authors

Bedelbaeva K, Snyder A, Gourevitch D, Clark L, Zhang XM, Leferovich J, Cheverud JM, Liet

#### Journal

Proc Natl Acad Sci U S A (2010) 107:5845-5850

#### Product Usage

Flow cytometric analysis of cell cycle progression in fibroblast-like cells from MRL mice.

#### ID:

PN68689

- [Communication via gap junctions underlies early functional and beneficial interactions b](#)  
[host.](#)

#### Authors

Jäderstad J, Jäderstad LM, Li J, Chintawar S, Salto C, Pandolfo M, Ourednik V, Teng YD, Sidi

#### Journal

Proc Natl Acad Sci U S A (2010) 107:5184-5189

#### Product Usage

Analysis of gap junction connectivity in cultured neural stem cells.

#### ID:

PN68236



- **Conventional apoptosis assays using propidium iodide generate a significant number of false assessment of cell death.**

Authors

Rieger AM, Hall BE, Luong le T, Schang LM, Barreda DR,

Journal

J Immunol Methods (2010) 358:81-92

Product Usage

Assessment of false positive results in annexin V/propidium iodide (PI) double labeling protocol for RNA

ID:

PN68179

- **High-content screening for biofilm assays.**

Authors

Peng F, Hoek EM, Damoiseaux R,

Journal

J Biomol Screen (2010) 15:748-754

Product Usage

High-content imaging assessment of biofilm formation and removal on engineered surfaces.

ID:

PN68697

- **Analysis of plasma membrane integrity by fluorescent detection of Tl(+) uptake.**

Authors

Bowman AM, Nesin OM, Pakhomova ON, Pakhomov AG,

Journal

J Membr Biol (2010) 236:15-26

Product Usage

Detection of membrane nanopores produced by ultrashort electric pulses.

ID:

PN68287

## Product Reviews (0)

## Help Center



### **Order Support**

Search frequently asked questions, register as a customer, change your account settings, track your orders, report a problem, or simply find out how to contact us.



### **Technical Support**

Tools and resources for the Support team. Find product guides, troubleshooting



People Who Viewed  
This Item Also Viewed  
[SYTOX® AADvanced™ Dead Cell Stain Kit](#)

Cat #

S10274

Size

1 kit

List Price

**325.50**

(USD)



[LIVE/DEAD® Fixable Red Dead Cell Stain Kit \\*for 488 nm excitation\\* \\*200 assays\\*](#)

Cat #

L23102

Size

1 kit

List Price

**243.00**

(USD)



[SYTO® 60 red fluorescent nucleic acid stain - 5 mM solution in DMSO](#)

Cat #

S11342

Size

250 µl

List Price

**189.00**

(USD)

## Product Categories

- [Nucleic Acid Cell Stains](#)

## Related Applications

- [Cell Structure](#)
- [Cell Tracing & Tracking](#)
- [Cell Viability, Proliferation & Function](#)
- [Cellular Imaging & Microscopy](#)
- [Chromatin Biology](#)

[View More](#)

## Related Products

- [Nuclear Probes](#)
- [Structural Probes](#)



## Fluorescent and Biotinylated Dextran—Section 14.5

### Page Contents

[Properties of Molecular Probes Dextran Conjugates](#)  
[Loading Cells with Dextran and Subsequent Tissue Processing](#)  
[Neuronal Tracing with Dextran](#)  
[Cell Lineage Tracing with Dextran](#)  
[Studying Intercellular Communication with Dextran](#)  
[Probing Membrane Permeability with Dextran](#)  
[Following Endocytosis with Dextran](#)  
[Tracing Fluid Transport with Dextran](#)  
[Ordering Information](#)

### Related Tables and Notes

[Techniques for loading molecules into the cytoplasm—Table 14.1](#)  
[Fluorescence characteristics of NeuroTrace fluorescent Nissl stains—Table 14.2](#)  
[Summary of Molecular Probes lipophilic carboxyfluorescein and aminostyryl tracers—Table 14.3](#)  
[Molecular Probes dextran conjugates—Table 14.4](#)  
[FluoSpheres microspheres for blood flow determination—Table 14.5](#)  
[FluoSpheres blood flow and color kits—Table 14.6](#)  
[FluoSpheres and TransFluoSpheres microspheres for tracing—Table 14.7](#)  
[Molecular Probes europium and platinum luminescent FluoSpheres microspheres—Table 14.8](#)  
[Anti-Lucifer Yellow Dye, Anti-Alexa Fluor 405/Cascade Blue Dye and Anti-Alexa Fluor 488 Dye Antibodies—Note 14.1](#)  
[Fluorescent Probes for Photoconversion of Diaminobenzidine Reagents—Note 14.2](#)  
[Assays of Volume Change, Membrane Fusion and Membrane Permeability—Note 14.3](#)

Share |

Dextran is a hydrophilic polysaccharide characterized by their moderate to high molecular weight, good water solubility and low toxicity. They are widely used as both anterograde and retrograde tracers in neurons [\(REF\)](#) and for many other applications. Dextran is biologically inert due to their uncommon poly-( $\alpha$ -D-1,6-glucose) linkages, which render them resistant to cleavage by most endogenous cellular glycosidases. They also usually have low immunogenicity.

We offer almost 100 fluorescent and biotinylated dextran conjugates in several molecular weight ranges. Because the source and molecular weight of the dextran, as well as the net charge, degree of substitution and nature of the dye may significantly affect the application, references citing the use of Molecular Probes dextran may not be directly applicable to dextran obtained from other sources and should be considered guides rather than definitive protocols. In most cases, Molecular Probes fluorescent dextran are much brighter and have higher negative charge than dextran available from other sources. Furthermore, we use rigorous methods for removing as much unconjugated dye as practical, and then assay our dextran conjugates by thin-layer chromatography to ensure the absence of low molecular weight contaminants.

### Properties of Molecular Probes Dextran Conjugates

#### A Wide Selection of Substituents

Molecular Probes dextran are conjugated to biotin or a wide variety of fluorophores, including seven of our Alexa Fluor dyes ([Molecular Probes dextran conjugates—Table 14.4](#)). In particular, we would like to highlight the dextran conjugates of Alexa Fluor 488, Oregon Green and Rhodamine Green dyes, which are significantly brighter and more photostable than most fluorescein dextran. Dextran-conjugated fluorescent indicators for calcium and magnesium ions ([Fluorescent  \$Ca^{2+}\$  Indicator Conjugates—Section 19.4](#)) and for pH ([pH Indicator Conjugates—Section 20.4](#)) are described with their corresponding ion indicators in other chapters.

#### Dextran Size

Molecular Probes dextran include those with nominal molecular weights (MW) of 3000; 10,000; 40,000; 70,000; 500,000; and 2,000,000 daltons ([Molecular Probes dextran conjugates—Table 14.4](#)). Because unlabeled dextran are polydisperse—and may become more so during the chemical processes required for their modification and purification—the actual molecular weights present in a particular sample may have a broad distribution. For example, our “3000 MW” dextran preparations contain polymers with molecular weights predominantly in the range of ~1500–3000 daltons, including the dye or other label.

#### Degree of Substitution of Molecular Probes Dextran

Dextran from other commercial sources usually have a degree of substitution of 0.2 or fewer dye molecules per dextran molecule for dextran in the 10,000 MW range. Molecular Probes dextran, however, typically contain 0.3–0.7 dyes per dextran in the 3000 MW range, 0.5–2 dyes per dextran in the 10,000 MW range, 2–4 dyes in the 40,000 MW range and 3–6 dyes in the 70,000 MW range. The actual degree of substitution is indicated on the product's label. If too many fluorophores are conjugated to the dextran molecule, quenching and undesired interactions with cellular components may occur. We have found our degree of substitution to be optimal for most applications, yielding dextran that are typically much more fluorescent than the labeled dextran available from other sources, thus permitting lower quantities to be used for intracellular tracing.


It has been reported that some commercially available fluorescein isothiocyanate (FITC) dextran yield spurious results in endocytosis studies because of the presence of free dye or metal contamination. [\(REF\)](#) To overcome this problem, we remove as much of the free dye as possible by a combination of precipitation, dialysis, gel filtration and other techniques. The fluorescent dextran is then assayed by thin-layer chromatography (TLC) to ensure that it is free of low molecular weight dyes. We prepare several unique products that have two or even three different labels, including the fluoro-ruby, mini-ruby and micro-ruby dextran, described below. Not all of the individual dextran molecules of these products are expected to have all the substituents, or to be equally fixable, particularly in conjugates of the lowest molecular weight dextran.






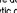
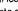
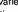
## Dextran Net Charge and Method of Substitution



The net charge on the dextran depends on the fluorophore and the method of preparing the conjugate. We prepare most of Molecular Probes dextran conjugates by reacting a water-soluble amino dextran (D1850, D1861, D1862, D3330, D7144) with the succinimidyl ester of the appropriate dye, rather than reacting a native dextran with isothiocyanate derivatives such as FITC. This method provides superior amine selectivity and yields an amide linkage, which is somewhat more stable than the corresponding thioureas formed from isothiocyanates. Except for the Rhodamine Green and Alexa Fluor 488 conjugates, once the dye has been added, the unreacted amines on the dextran are capped to yield a neutral or anionic dextran. In the case of the Rhodamine Green and Alexa Fluor 488 dextran conjugates, the unreacted amines on the dextran are not capped after dye conjugation. Thus, these dextran conjugates may be neutral, anionic or cationic. The Alexa Fluor, Cascade Blue, lucifer yellow, fluorescein and Oregon Green dextran conjugates are intrinsically anionic, whereas most of the dextran labeled with the zwitterionic rhodamine B, tetramethylrhodamine and Texas Red dyes are essentially neutral. To produce more highly anionic dextran, we have developed a proprietary procedure for adding negatively charged groups to the dextran carriers; these products are designated "polyanionic" dextran.



## Dextran Fixability

Some applications require that the dextran tracer be treated with formaldehyde or glutaraldehyde for subsequent analysis.  For these applications, we offer "lysine-fixable" versions of most of our dextran conjugates of fluorophores or biotin. These dextran have covalently bound lysine residues that permit dextran tracers to be conjugated to surrounding biomolecules by aldehyde-mediated fixation for subsequent detection by immunohistochemical and ultrastructural techniques. We have also shown that all of our 10,000 MW Alexa Fluor dextran conjugates can be fixed with aldehyde-based fixatives; however, due to their smaller size, our Alexa Fluor 3000 MW dextran conjugates most likely will not survive fixation procedures.





## Loading Cells with Dextran and Subsequent Tissue Processing

 Unless taken up by an endocytic process, dextran conjugates are membrane impermeant and usually must be loaded by relatively invasive techniques.  Techniques for loading molecules into the cytoplasm—Table 14.1. As with the lipophilic tracers in Tracers for Membrane Labeling—Section 14.4, crystals of the dextran conjugates have been successfully loaded by simply placing them directly on some kinds of samples.  We have found the Influx pinocytotic cell-loading reagent (I14492). Chelators, Calibration Buffers, Ionophores and Cell-loading Reagents—Section 12.8 to be useful for loading dextran into a variety of adherent and nonadherent cells.  Sterile filtration of dextran solutions before use with live cells is highly recommended.  Biotin and biotinylated biomolecules with molecular weights up to >100,000 daltons are taken up by some plant cells through an endocytic pathway. 

Our lysine-fixable dextran and 10,000 MW Alexa Fluor dextran can be fixed in place with formaldehyde or glutaraldehyde, permitting subsequent tissue processing, such as sectioning. A protocol has been published for embedding tissues in plastic for high-resolution characterization of neurons filled with lysine-fixable fluorescent dextran.  Fixation of biotinylated or fluorescent dextran also permits the use of fluorescent- or enzyme-labeled conjugates of avidin and streptavidin (Avidin, Streptavidin, NeutrAvidin and CaptAvidin Biotin-Binding Proteins and Affinity Matrices—Section 7.6, Molecular Probes avidin, streptavidin, NeutrAvidin and CaptAvidin conjugates—Table 7.23) or of anti-dye antibodies (Anti-Dye and Anti-Hapten Antibodies—Section 7.4, Anti-fluorophore and anti-hapten antibodies—Table 7.18), respectively. These techniques can amplify the signal, which is important for detecting fine structure in sections or for changing the detection mode.  We provide antibodies to the Alexa Fluor 488, Alexa Fluor 405/Cascade Blue, lucifer yellow, fluorescein, BODIPY FL, tetramethylrhodamine and Texas Red fluorophores and to the 2,4-dinitrophenyl (DNP) and nitrotyrosine haptens (Anti-Dye and Anti-Hapten Antibodies—Section 7.4).


Photocrosslinking of neurons labeled with lysine-fixable fluorescent dextran in the presence of diaminebenzidine (DAB) using the Diaminebenzidine (DAB) Histochemistry Kits (Secondary Immunoreagents—Section 7.2, Avidin, Streptavidin, NeutrAvidin and CaptAvidin Biotin-Binding Proteins and Affinity Matrices—Section 7.6) can be used to produce electron-dense products for electron microscopy.  (Fluorescent Probes for Photocrosslinking of Diaminebenzidine Reagents—Note 14.2). Electron-dense products can also be generated from peroxidase or colloidal gold conjugates of avidin, streptavidin or anti-dye antibodies.  NANOGOLD and Alexa Fluor FluoroNanogold conjugates of secondary antibodies (Secondary Immunoreagents—Section 7.2) and streptavidin (Avidin, Streptavidin, NeutrAvidin and CaptAvidin Biotin-Binding Proteins and Affinity Matrices—Section 7.6) can be utilized to allow detection of labeled dextran in fixed-cell preparations by light microscopy or, following silver enhancement with the LI Silver Enhancement Kit (L24919, Secondary Immunoreagents—Section 7.2), by electron microscopy.






## Neuronal Tracing with Dextran


Fluorescent and biotinylated dextran are routinely employed to trace neuronal projections. Dextran can function efficiently as anterograde or retrograde tracers,  depending on the study method and tissue type used. Active transport of dextran occurs only in live, not fixed tissue.  Comparative studies of rhodamine isothiocyanate, rhodamine B dextran (D1824) and lysinated tetramethylrhodamine dextran (fluoro-ruby, D1817) have shown that the dextran conjugates produce less diffusion at injection sites and more permanent labeling than do the corresponding free dyes.  Dextran conjugates with molecular weights up to 70,000 daltons have been employed as neuronal tracers in a wide variety of species. The availability of fluorescent dextran conjugates with different sizes and charges permitted the analysis of direction and rate of axonal transport in the squid giant axon. 

## Multilabeled Dextran

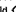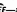

Molecular Probes fixable dextran, most of which are lysinated dextran (see the products marked by a single dagger (†) in Molecular Probes dextran conjugates—Table 14.4), are generally preferred for neuronal tracing because they may transport more effectively and can be fixed in place with aldehydes after labeling. We prepare a number of multilabeled dextran that are fixable, including some that have acquired the distinction of unique names in various publications:

**Fluoro-ruby** —a red-orange-fluorescent, aldehyde-fixable 10,000 MW dextran labeled with both tetramethylrhodamine and lysine (D1817). 3000 MW, 70,000 MW and 2,000,000 MW versions of fluoro-ruby are also available (D3308, D1818, D7139).

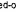**Fluoro-emerald** —a green-fluorescent, aldehyde-fixable 10,000 MW dextran labeled with both fluorescein and lysine (D1820;   ). This labeled dextran is also available in molecular weights from 3000 daltons up to 2,000,000 daltons  (D3306, D1845, D1822, D7136, D7137).

**Micro-ruby** (D7182) and **mini-ruby**  (D3312)—red-orange-fluorescent, aldehyde-fixable 3000 MW and 10,000 MW dextran simultaneously labeled with tetramethylrhodamine, biotin and lysine.

**Micro-emerald** (D7156) and **mini-emerald** (D7178)—green-fluorescent, aldehyde-fixable dextran simultaneously labeled with fluorescein, biotin and lysine.

**Biotinylated dextran amine (BDA)** —nonfluorescent, aldehyde-fixable dextran simultaneously labeled with biotin and lysine and available in several molecular weights (D1856, D1957, D7135, D7142). A useful review has been published on the BDA products. 



Fluoro-ruby and fluoro-emerald (**REF**) have been extensively employed for retrograde and anterograde neuronal tracing, **REF** transplantation **REF** and cell-lineage tracing. **REF** Both products can be used to photoconvert DAB into an insoluble, electron-dense reaction product. **REF** Like fluoro-ruby and fluoro-emerald, micro-ruby and mini-ruby are brightly fluorescent, making it easy to visualize the electrode during the injection process. **REF** **Tracers for Membrane Labeling—Section 14.4**) or other lipophilic probes in **Tracers for Membrane Labeling—Section 14.4** can be used to mark the sites of microinjection. **REF** In addition, because these dextrans include a covalently linked biotin, filled cells can be probed with standard enzyme-labeled avidin or streptavidin conjugates or with NANOGOLD and Alexa Fluor FluoroNanogold streptavidin (**Avidin, Streptavidin, NeutrAvidin and CapAvidin Biotin-Binding Proteins and Affinity Matrices—Section 7.6**) to produce a permanent record of the experiment. **REF** Mini-ruby has proven useful for intracellular filling in fixed brain slices **REF** and has been reported to produce staining comparable to that achieved with lucifer yellow CH **REF** (**L453, L692, L1177, L12926, Polar Tracers—Section 14.3**). Moreover, the use of mini-ruby in conjunction with standard peroxidase-mediated avidin-biotin methods does not cause co-conversion of lipofuscin granules found in adult human brain, a common problem during photoconversion of lucifer yellow CH. **REF** The lysine-fixable micro-emerald and mini-emerald dextrans (triply labeled with fluorescein, biotin and lysine) provide a contrasting color that is better excited by the argon-ion laser of confocal laser-scanning microscopes; they have uses similar to micro-ruby and mini-ruby, respectively.

### 3000 MW Dextrans

The nominally 3000 MW dextrans offer several advantages over higher molecular weight dextrans, including faster axonal diffusion and greater access to peripheral cell processes **REF** (**REF**). Our "3000 MW" dextran preparations contain polymers with molecular weight predominantly in the range of ~1500–3000 daltons, including the dye or other label. Our selection of 3000 MW dextrans includes Alexa Fluor, fluorescein, Rhodamine Green, tetramethylrhodamine, Texas Red and biotin conjugates. We also offer lysine-fixable 3000 MW dextrans that are simultaneously labeled with both fluorescein and biotin (micro-emerald, **D7155**) or tetramethylrhodamine and biotin (micro-ruby, **D7162**).

The 3000 MW fluorescein dextran and tetramethylrhodamine dextran (**D3306, D3308; REF, REF**) have been observed to readily undergo both anterograde and retrograde movement in live cells. **REF** These 3000 MW dextrans appear to passively diffuse within the neuronal process, as their intracellular transport is not effectively inhibited by colchicine or nocodazole, both of which disrupt active transport by depolymerizing microtubules. **REF** Moreover, these small dextrans diffuse at rates equivalent to those of smaller tracers such as sulfohodamine 101 and biocytin (~2 millimeters/hour at 22°C) and about twice as fast as 10,000 MW dextrans. The relatively low molecular weight of the dextrans may result in transport of some labeled probes through gap junctions (see below). The signal from tetramethylrhodamine-conjugated dextrans can be detected in the fine dendritic configuration of cortical projection neurons using anti-tetramethylrhodamine antibodies (**A6387, Anti-Dye and Anti-Hapten Antibodies—Section 7.4**) and peroxidase-anti-peroxidase complex staining. **REF**

### NeuroTrace BDA-10,000 Neuronal Tracer Kit

The NeuroTrace BDA-10,000 Neuronal Tracer Kit (**N7167**) contains convenient amounts of each of the components required for neuroanatomical tracing using BDA methods, **REF** including:

- Lysine-fixable, biotinylated 10,000 MW dextran amine (BDA-10,000)
- Horse radish peroxidase avidin (HRP avidin)
- 3,3'-Diaminobenzidine (DAB)
- Rigorously tested protocols for fast and simple tracing experiments (**NeuroTrace BDA-10,000 Neuronal Tracer Kit**)

The neuronal tracer BDA-10,000 is transported over long distances and fills fine processes bidirectionally, including boutons in the anterograde direction and dendritic structures in the retrograde direction. **REF** Two days to two weeks after BDA-10,000 is injected into the desired region of the brain, the brain tissue can be fixed and sectioned. BDA-10,000 can also be applied to cut nerves and allowed to transport. Following incubation with HRP avidin and then DAB, the electron-dense DAB reaction product can be viewed by either light or electron microscopy **REF** (**REF**). The NeuroTrace BDA-10,000 labeling method can be readily combined with other anterograde or retrograde labeling methods or with immunohistochemical techniques. BDA-10,000 is available as a separate product (**D1958**), as are BDA derivatives with other molecular weights—BDA-3000 **REF** (**D7135**), BDA-7000 **REF** (**D1957**) and BDA-500,000 **REF** (**D7142**). A detailed protocol that utilizes Molecular Probes BDA-10,000 probe, HRP streptavidin (**S911, Avidin, Streptavidin, NeutrAvidin and CapAvidin Biotin-Binding Proteins and Affinity Matrices—Section 7.6**) and tetramethylbenzidine to anterogradely label fine processes in neurons has been published. **REF**

### Cell Lineage Tracing with Dextrans

Fluorescent dextrans—particularly the fluorescein and rhodamine conjugates—have been used extensively for tracing cell lineage. **REF** Our Alexa Fluor 647 and Alexa Fluor 680 dextrans (**D22214, D34680, D34681**) provide longer-wavelength detection options for specimens with high levels of autofluorescence or low transparency. In this technique, the dextran is microinjected into a single cell of the developing embryo, and the fate of that cell and its daughters can be followed *in vivo* (**REF**). The lysine-fixable tetramethylrhodamine and Texas Red dextran conjugates (**Molecular Probes dextran conjugates—Table 14.4**) are most frequently cited for lineage tracing studies. As a second color, particularly in combination with the Texas Red dextrans, researchers have most often used Molecular Probes lysine-fixable fluorescein dextrans **REF** (e.g., **D3306, D1820, D1822**). Although these fixable conjugates can be employed with long-term preservation of the tissue, some researchers prefer to co-inject a fluorescent, nonsynthesized dextran along with a nonfluorescent, lysine-fixable biotin dextran (BDA). The nonfluorescent BDA can then be fixed in place with aldehyde-based fixatives and probed with any of our fluorescent or enzyme-labeled streptavidin and avidin conjugates described in **Avidin, Streptavidin, NeutrAvidin and CapAvidin Biotin-Binding Proteins and Affinity Matrices—Section 7.6** (**Molecular Probes avidin, streptavidin, NeutrAvidin and CapAvidin conjugates—Table 7.23**).

Our 500,000 and 2,000,000 MW fluorescent dextrans (**Molecular Probes dextran conjugates—Table 14.4**) may be particularly useful for lineage tracing at early stages of development, although these biopolymers have lower water solubility and a greater tendency to precipitate or clog microinjection needles than our lower molecular weight dextrans. Some studies suggest that lower molecular weight dextrans may leak from blastomeres, complicating analysis. Injection of 2,000,000 MW fluorescein dextran and 2,000,000 MW Texas Red dextran into separate cells of the two-cell stage zebrafish embryo allowed the construction of a fate map. **REF** The 500,000 MW and 2,000,000 MW dextrans are labeled with fluorescein, tetramethylrhodamine or Texas Red dyes or with biotin, and all contain aldehyde-fixable lysine groups. The nonfluorescent 500,000 MW aminodextran (**D7144**) can be conjugated with the researcher's choice of amine-reactive reagents. **REF**

### Studying Intercellular Communication with Dextrans

The size of dextrans can be exploited to study connectivity between cells. **REF** Examples include studies of the passage of 3000 MW dextrans through plasmodesmata **REF** and modulation of gap junctional communication by transforming growth factor- $\beta$  and forskolin. **REF** However, the dispersion of molecular weights in our "3000 MW" dextran preparations, which contain polymers with total molecular weights predominantly in the range of ~1500–3000 daltons but may also contain molecules <1500 daltons, may complicate such analyses.



An important experimental approach to identifying cells that form gap junctions makes use of simultaneous introduction of the polar tracer lucifer yellow CH (~450 daltons) and a tetramethylrhodamine 10,000 MW dextran. Because low molecular weight tracers like lucifer yellow CH (L453, L12926; [Polar Tracers—Section 14.3](#)) pass through gap junctions and dextrans do not, the initially labeled cell exhibits red fluorescence, whereas cells connected through gap junctions have yellow fluorescence ([REF](#) (Figure 14.79)). This technique has been used to follow the loss of intercellular communication in adenocarcinoma cells, ([REF](#)) to show the re-establishment of communication during wound healing in *Drosophila* ([REF](#)) and to investigate intercellular communication at different stages in *Xenopus* embryos. ([REF](#)) Simultaneous loading of cells with two (or more) dextrans that differ in both their molecular weight and in the dye's fluorescence properties has been used to assess subcellular heterogeneities in the submicroscopic structure of cytoplasm. ([REF](#))

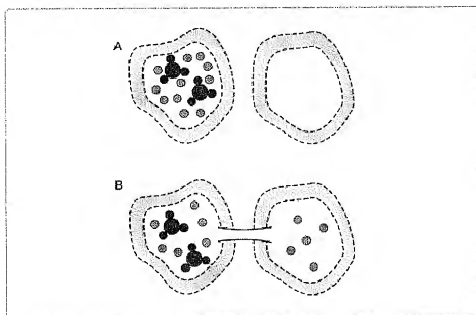


Figure 14.79 Dual-tracer technique for identifying gap junction-coupled cells. A) Cells are labeled with a mixture of a small polar tracer such as lucifer yellow CH (green circles) and a relatively large tetramethylrhodamine-labeled dextran (red circles). B) Adjoining gap junction-coupled cells are accessible to the low molecular weight tracer whereas the much larger dextran conjugate is excluded. Coupled cells with single-color lucifer yellow CH labeling are readily distinguished from initially labeled cells with dual fluorescence.

### Probing Membrane Permeability with Dextrans

Labeled dextrans are often used to investigate vascular permeability and blood-brain barrier integrity. ([REF](#)) Fluorescein dextrans with molecular weights ranging from 4000 to 150,000 daltons have been used to determine the effect of electroporation variables—pulse size, shape and duration—on plasma-membrane pore size in chloplast, ([REF](#)) red blood cells ([REF](#)) and fibroblasts. ([REF](#)) Fluorescence recovery after photobleaching (FRAP) techniques have been used to monitor nucleocytoplasmic transport of fluorescent dextrans of various molecular weights, allowing the determination of the size-exclusion limit of the nuclear pore membrane, ([REF](#)) as well as to study the effect of epidermal growth factor and insulin on the nuclear membrane and on nucleocytoplasmic transport. ([REF](#))

Microinjected 3000 MW fluorescent dextrans concentrate in interphase nuclei of *Drosophila* embryos, whereas 40,000 MW dextrans remain in the cytoplasm and enter the nucleus only after breakdown of the nuclear envelope during prophase. This size-exclusion phenomenon was used to follow the cyclical breakdown and reformation of the nuclear envelope during successive cell divisions. ([REF](#)) Similarly, our 10,000 MW Calcium Green dextran conjugate (C3713; [Fluorescent Ca<sup>2+</sup> Indicator Conjugates—Section 19.4](#)) was shown to diffuse across the nuclear membrane of isolated nuclei from *Xenopus laevis* oocytes, but the 70,000 MW and 500,000 MW conjugates could not. ([REF](#)) Significantly, depletion of nuclear Ca<sup>2+</sup> stores by inositol 1,4,5-trisphosphate (Ins 1,4,5-P<sub>3</sub>, I3716; [Calcium Regulation—Section 17.2](#)) or by calcium chelators ([Chelators, Calibration Buffers, Ionophores and Cell-Loading Reagents—Section 19.8](#)) blocked nuclear uptake of the 10,000 MW Calcium Green dextran conjugate but not entry of lucifer yellow CH. Our 3000 MW Calcium Green dextran conjugate (C8785) is actively transported in adult nerve fibers over a significant distance and is retained in presynaptic terminals in a form that allows monitoring of presynaptic Ca<sup>2+</sup> levels. ([REF](#))

### Following Endocytosis with Dextrans

#### Fluorescent Dextrans

Tracing internalization of extracellularly introduced fluorescent dextrans is a standard method for analyzing fluid-phase endocytosis. ([REF](#)) We offer dextrans with nominal molecular weights ranging from 3000 to 2,000,000 daltons, many of which can also be used as pinocytosis or phagocytosis markers ([Molecular Probes dextran conjugates—Table 14.4](#)). Discrimination of internalized fluorescent dextrans from dextrans in the growth medium is facilitated by use of reagents that quench the fluorescence of the external probe. For example, most of our anti-fluorophore antibodies ([Anti-Dye and Anti-Hapten Antibodies—Section 7.4](#), [Anti-fluorophore and anti-hapten antibodies—Table 7.19](#)) strongly quench the fluorescence of the corresponding dyes. ([REF](#))

Negative staining produced by fluorescent dextrans that have been intracellularly infused via a patch pipette is indicative of nonendocytic vacuoles in live pancreatic acinar cells; extracellular addition of a second, color-contrasting dextran then allows discrimination of endocytic and nonendocytic vacuoles. ([REF](#)) An *in vitro* assay for homotypic fusion of early endosomes has been described in which two cell populations are labeled with Alexa Fluor 488 and Alexa Fluor 594 10,000 MW dextrans (D22910, D22914) by fluid phase uptake, followed by subcellular fractionation and analysis of endosomal fluorescence colocalization. ([REF](#)) Intracellular fusion of endosomes has also been monitored with a BODIPY FL avidin conjugate by following the fluorescence enhancement that occurs when it complexes with a biotinylated dextran. ([REF](#))



## pH Indicator Dextrans

Some of the dyes we use to prepare Molecular Probes dextran conjugates exhibit fluorescence that is sensitive to the pH of the medium; pH indicator dextrans and their optical responses are described in detail in [pH Indicator Conjugates—Section 20.4](#). Consequently, internalization of labeled dextran into acidic organelles of cells can often be tracked by measuring changes in the fluorescence of the dye. The fluorescein dextrans ( $pK_a \sim 6.4$ ) are frequently used to investigate endocytic acidification. [REF](#) Fluorescence of fluorescein-labeled dextrans is strongly quenched upon acidification; however, fluorescein's lack of a spectral shift in acidic solution makes it difficult to discriminate between an internalized probe that is quenched and residual fluorescence of the external medium. Dextran conjugates that either shift their emission spectra in acidic environments, such as the SNARF dextrans ([pH Indicator Conjugates—Section 20.4](#)), or undergo significant shifts of their excitation spectra, such as BCECF and Oregon Green dextrans ([pH Indicator Conjugates—Section 20.4](#)), provide alternatives to fluorescein. The Oregon Green 488 and Oregon Green 514 dextrans exhibit a  $pK_a$  of approximately 4.7, facilitating measurements in acidic environments. [REF](#) In addition to these pH indicator dextrans, we prepare a dextran that is double-labeled with fluorescein and tetramethylrhodamine ([D1851: pH Indicator Conjugates—Section 20.4](#)), which has been used as a ratiometric indicator (Figure 20.28) to measure endosomal acidification in Hep G2 cells [REF](#) and murine alveolar macrophages. [REF](#)

In contrast to fluorescein and Oregon Green 488 dextrans, pHrodo 10,000 MW dextran ([P10361](#)) exhibits increasing fluorescence in response to acidification [REF](#) (Figure 16.47). The minimal fluorescent signal from pHrodo dextran at neutral pH prevents the detection of noninternalized and nonspecifically bound conjugates and eliminates the need for quenching reagents and extra wash steps, thus providing a simple fluorescent assay for endocytic activity. pHrodo dextran's excitation and emission maxima of 560 and 585 nm, respectively, facilitate multiplexing with other fluorophores including blue-, green- and far-red-fluorescent probes. Although pHrodo dextran is optimally excited at approximately 560 nm, it is also readily excited by the 488 nm spectral line of the argon-ion laser found on flow cytometers, confocal microscopes and imaging microplate readers (Figure 16.48).

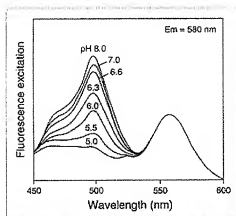


Figure 20.28 The excitation spectra of double-labeled fluorescein-tetramethylrhodamine dextran ([D1851](#)), which contains pH-dependent (fluorescein) and pH-independent (tetramethylrhodamine) dyes.

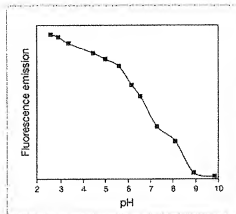


Figure 16.47 The pH response profile of pHrodo dextran ([P10361](#)) monitored at excitation/emission wavelengths of 545/590 nm in a fluorescence microplate reader. Citrate, MOPS and borate buffers were used to span the pH range from 2.5 to 10.

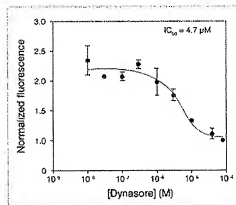


Figure 16.48 Tracking endocytosis inhibition with pHrodo dextran conjugates. HeLa cells were plated in 96-well format and treated with dynasore for 3 hours at 37°C prior to the pHrodo endocytosis assay. Next, 40  $\mu\text{g}/\text{mL}$  of pHrodo 10,000 MW dextran ([P10361](#)) was incubated for 30 minutes at 37°C, and cells were then stained with HCS NuclearMask Blue Stain ([H10325](#)) for 10 minutes to reveal total cell number and demarcation for image analysis. Images were acquired on the BD Pathway 855 High-Content Biomager (BD Biosciences).



Fluorescent dextrans are important tools for studying the hydrodynamic properties of the cytoplasmic matrix. The intracellular mobility of these fluorescent tracers can be investigated using fluorescence recovery after photobleaching (FRAP) techniques. We offer a range of dextran sizes, thus providing a variety of hydrodynamic radii for investigating both the nature of the cytoplasmic matrix and the permeability of the surrounding membrane. Because of their solubility and biocompatibility, fluorescent dextrans have been used to monitor *in vivo* tissue permeability and flow in the uveoscleral tract, *capillaries* and proximal tubules, as well as diffusion of high molecular weight substances in the brain's extracellular environment.

## Ordering Information

SKU	Product Catalog	Size	Price	Quantity
D-1816	dextran, tetramethylrhodamine, 10,000 MW, neutral	25 mg	140.00 USD	<input type="text"/>
D-1817	dextran, tetramethylrhodamine, 10,000 MW, lysine fixable (fluoro-ruby)	25 mg	189.00 USD	<input type="text"/>
D-1818	dextran, tetramethylrhodamine, 70,000 MW, lysine fixable	25 mg	189.00 USD	<input type="text"/>
D-1819	dextran, tetramethylrhodamine, 70,000 MW, neutral	25 mg	140.00 USD	<input type="text"/>
D-1820	dextran, fluorescein, 10,000 MW, anionic, lysine fixable (fluoro-emerald)	25 mg	189.00 USD	<input type="text"/>
D-1821	dextran, fluorescein, 10,000 MW, anionic	25 mg	140.00 USD	<input type="text"/>
D-1822	dextran, fluorescein, 70,000 MW, anionic, lysine fixable	25 mg	189.00 USD	<input type="text"/>
D-1823	dextran, fluorescein, 70,000 MW, anionic	25 mg	140.00 USD	<input type="text"/>
D-1824	dextran, rhodamine B, 10,000 MW, neutral	25 mg	140.00 USD	<input type="text"/>
D-1825	dextran, lucifer yellow, 10,000 MW, anionic, lysine fixable	25 mg	140.00 USD	<input type="text"/>
D-1826	dextran, Texas Red®, 10,000 MW, neutral	25 mg	142.00 USD	<input type="text"/>
D-1828	dextran, Texas Red®, 40,000 MW, neutral	25 mg	140.00 USD	<input type="text"/>
D-1830	dextran, Texas Red®, 70,000 MW, neutral	25 mg	140.00 USD	<input type="text"/>
D-1841	dextran, rhodamine B, 70,000 MW, neutral	25 mg	140.00 USD	<input type="text"/>
D-1842	dextran, tetramethylrhodamine, 40,000 MW, neutral	25 mg	140.00 USD	<input type="text"/>
D-1844	dextran, fluorescein, 40,000 MW, anionic	25 mg	140.00 USD	<input type="text"/>
D-1845	dextran, fluorescein, 40,000 MW, anionic, lysine fixable	25 mg	189.00 USD	<input type="text"/>
D-1860	dextran, amino, 10,000 MW	1 g	131.00 USD	<input type="text"/>
D-1861	dextran, amino, 40,000 MW	1 g	131.00 USD	<input type="text"/>
D-1862	dextran, amino, 70,000 MW	1 g	131.00 USD	<input type="text"/>
D-1863	dextran, Texas Red®, 10,000 MW, lysine fixable	25 mg	189.00 USD	<input type="text"/>
D-1864	dextran, Texas Red®, 70,000 MW, lysine fixable	25 mg	189.00 USD	<input type="text"/>
D-1868	dextran, tetramethylrhodamine, 10,000 MW, anionic, fixable	25 mg	140.00 USD	<input type="text"/>
D-1958	dextran, biotin, 10,000 MW, lysine fixable (BDA-10,000)	25 mg	189.00 USD	<input type="text"/>
D-1957	dextran, biotin, 70,000 MW, lysine fixable (BDA-70,000)	25 mg	189.00 USD	<input type="text"/>
D-1976	dextran, Cascade Blue®, 10,000 MW, anionic, lysine fixable	25 mg	189.00 USD	<input type="text"/>
D-22910	dextran, Alexa Fluor® 488, 10,000 MW, anionic, fixable	5 mg	258.00 USD	<input type="text"/>
D-22911	dextran, Alexa Fluor® 546, 10,000 MW, anionic, fixable	5 mg	258.00 USD	<input type="text"/>
D-22912	dextran, Alexa Fluor® 568, 10,000 MW, anionic, fixable	5 mg	258.00 USD	<input type="text"/>
D-22913	dextran, Alexa Fluor® 594, 10,000 MW, anionic, fixable	5 mg	258.00 USD	<input type="text"/>
D-22914	dextran, Alexa Fluor® 647, 10,000 MW, anionic, fixable	2 mg	140.00 USD	<input type="text"/>
D-3305	dextran, fluorescein, 3000 MW, anionic	10 mg	189.00 USD	<input type="text"/>
D-3306	dextran, fluorescein, 3000 MW, anionic, lysine fixable	10 mg	214.00 USD	<input type="text"/>



D-3307	dextran, tetramethylrhodamine, 3000 MW, anionic	10 mg	189.00 USD	<input type="text"/>
D-3308	dextran, tetramethylrhodamine, 3000 MW, anionic, lysine fixable	10 mg	214.00 USD	<input type="text"/>
D-3312	dextran, tetramethylrhodamine and biotin, 10,000 MW, lysine fixable (mini-ruby)	10 mg	258.00 USD	<input type="text"/>
D-3328	dextran, Texas Red®, 3000 MW, lysine fixable	10 mg	214.00 USD	<input type="text"/>
D-3329	dextran, Texas Red®, 3000 MW, neutral	10 mg	189.00 USD	<input type="text"/>
D-3330	dextran, amino, 3000 MW	100 mg	95.50 USD	<input type="text"/>
D34679	dextran, Alexa Fluor® 555, 10,000 MW, anionic, fixable	5 mg	274.00 USD	<input type="text"/>
D34680	dextran, Alexa Fluor® 680, 10,000 MW, anionic, fixable	5 mg	274.00 USD	<input type="text"/>
D34681	dextran, Alexa Fluor® 680, 3,000 MW, anionic	2 mg	160.00 USD	<input type="text"/>
D34682	dextran, Alexa Fluor® 488, 3,000 MW, anionic	2 mg	160.00 USD	<input type="text"/>
D-7132	dextran, Cascade Blue®, 3000 MW, anionic, lysine fixable	10 mg	214.00 USD	<input type="text"/>
D-7135	dextran, biotin, 3000 MW, lysine fixable (BDA-3000)	10 mg	214.00 USD	<input type="text"/>
D-7138	dextran, fluorescein, 500,000 MW, anionic, lysine fixable	10 mg	189.00 USD	<input type="text"/>
D-7137	dextran, fluorescein, 2,000,000 MW, anionic, lysine fixable	10 mg	189.00 USD	<input type="text"/>
D-7138	dextran, tetramethylrhodamine, 2,000,000 MW, lysine fixable	10 mg	189.00 USD	<input type="text"/>
D-7142	dextran, biotin, 500,000 MW, lysine fixable (BDA-500,000)	10 mg	189.00 USD	<input type="text"/>
D-7144	dextran, amino, 500,000 MW	100 mg	95.50 USD	<input type="text"/>
D-7153	dextran, Rhodamine Green™, 10,000 MW, lysine fixable	10 mg	258.00 USD	<input type="text"/>
D-7156	dextran, fluorescein and biotin, 3000 MW, anionic, lysine fixable (micro-emerald)	5 mg	258.00 USD	<input type="text"/>
D-7162	dextran, tetramethylrhodamine and biotin, 3000 MW, lysine fixable (micro-ruby)	5 mg	258.00 USD	<input type="text"/>
D-7163	dextran, Rhodamine Green™, 3000 MW	5 mg	325.00 USD	<input type="text"/>
D-7168	dextran, BODIPY® FL, 10,000 MW, fixable	5 mg	189.00 USD	<input type="text"/>
D-7170	dextran, Oregon Green® 488, 10,000 MW, anionic	5 mg	189.00 USD	<input type="text"/>
D-7171	dextran, Oregon Green® 488, 10,000 MW, anionic, lysine fixable	5 mg	214.00 USD	<input type="text"/>
D-7172	dextran, Oregon Green® 488, 70,000 MW, anionic	5 mg	189.00 USD	<input type="text"/>
D-7173	dextran, Oregon Green® 488, 70,000 MW, anionic, lysine fixable	5 mg	214.00 USD	<input type="text"/>
D-7176	dextran, Oregon Green® 514, 70,000 MW, anionic	5 mg	189.00 USD	<input type="text"/>
D-7178	dextran, fluorescein and biotin, 10,000 MW, anionic, lysine fixable (mini-emerald)	10 mg	257.00 USD	<input type="text"/>
N-7167	NeuroTrace® BDA-10,000 Neuronal Tracer Kit	1 kit	259.00 USD	<input type="text"/>
P10361	dextran, phalloidin™, 10,000 MW for endocytosis	0.5 mg	265.00 USD	<input type="text"/>
V-22915	Vybrant® Cell Lineage Tracing Kit	1 kit	141.00 USD	<input type="text"/>

Add to Cart



## Ordering Information

SKU	Product Catalog	Size	Price
D-1816	<a href="#">dextran, tetramethylrhodamine, 10,000 MW, neutral</a>	25 mg	140.00 USD
D-1817	<a href="#">dextran, tetramethylrhodamine, 10,000 MW, lysine fixable (fluoro-ruby)</a>	25 mg	189.00 USD
D-1818	<a href="#">dextran, tetramethylrhodamine, 70,000 MW, lysine fixable</a>	25 mg	189.00 USD
D-1819	<a href="#">dextran, tetramethylrhodamine, 70,000 MW, neutral</a>	25 mg	140.00 USD
D-1820	<a href="#">dextran, fluorescein, 10,000 MW, anionic, lysine fixable (fluoro-emerald)</a>	25 mg	189.00 USD
D-1821	<a href="#">dextran, fluorescein, 10,000 MW, anionic</a>	25 mg	140.00 USD
D-1822	<a href="#">dextran, fluorescein, 70,000 MW, anionic, lysine fixable</a>	25 mg	189.00 USD
D-1823	<a href="#">dextran, fluorescein, 70,000 MW, anionic</a>	25 mg	140.00 USD
D-1824	<a href="#">dextran, rhodamine B, 10,000 MW, neutral</a>	25 mg	140.00 USD
D-1825	<a href="#">dextran, lucifer yellow, 10,000 MW, anionic, lysine fixable</a>	25 mg	140.00 USD
D-1828	<a href="#">dextran, Texas Red®, 10,000 MW, neutral</a>	25 mg	142.00 USD
D-1829	<a href="#">dextran, Texas Red®, 40,000 MW, neutral</a>	25 mg	140.00 USD
D-1830	<a href="#">dextran, Texas Red®, 70,000 MW, neutral</a>	25 mg	140.00 USD
D-1841	<a href="#">dextran, rhodamine B, 70,000 MW, neutral</a>	25 mg	140.00 USD
D-1842	<a href="#">dextran, tetramethylrhodamine, 40,000 MW, neutral</a>	25 mg	140.00 USD
D-1844	<a href="#">dextran, fluorescein, 40,000 MW, anionic</a>	25 mg	140.00 USD
D-1845	<a href="#">dextran, fluorescein, 40,000 MW, anionic, lysine fixable</a>	25 mg	189.00 USD
D-1860	<a href="#">dextran, amino, 10,000 MW</a>	1 g	131.00 USD
D-1861	<a href="#">dextran, amino, 40,000 MW</a>	1 g	131.00 USD
D-1862	<a href="#">dextran, amino, 70,000 MW</a>	1 g	131.00 USD
D-1863	<a href="#">dextran, Texas Red®, 10,000 MW, lysine fixable</a>	25 mg	189.00 USD
D-1864	<a href="#">dextran, Texas Red®, 70,000 MW, lysine fixable</a>	25 mg	189.00 USD
D-1868	<a href="#">dextran, tetramethylrhodamine, 10,000 MW, anionic, fixable</a>	25 mg	140.00 USD
D-1956	<a href="#">dextran, biotin, 10,000 MW, lysine fixable (BDA-10,000)</a>	25 mg	189.00 USD
D-1957	<a href="#">dextran, biotin, 70,000 MW, lysine fixable (BDA-70,000)</a>	25 mg	189.00 USD
D-1976	<a href="#">dextran, Cascade Blue®, 10,000 MW, anionic, lysine fixable</a>	25 mg	189.00 USD
D-22910	<a href="#">dextran, Alexa Fluor® 488, 10,000 MW, anionic, fixable</a>	5 mg	258.00 USD
D-22911	<a href="#">dextran, Alexa Fluor® 546, 10,000 MW, anionic, fixable</a>	5 mg	258.00 USD
D-22912	<a href="#">dextran, Alexa Fluor® 568, 10,000 MW, anionic, fixable</a>	5 mg	258.00 USD
D-22913	<a href="#">dextran, Alexa Fluor® 594, 10,000 MW, anionic, fixable</a>	5 mg	258.00 USD
D-22914	<a href="#">dextran, Alexa Fluor® 647, 10,000 MW, anionic, fixable</a>	2 mg	140.00 USD
D-3305	<a href="#">dextran, fluorescein, 3000 MW, anionic</a>	10 mg	189.00 USD



<u>D-3306</u>	<u>dextran, fluorescein, 3000 MW, anionic, lysine fixable</u>	10 mg	214.00 USD
<u>D-3307</u>	<u>dextran, tetramethylrhodamine, 3000 MW, anionic</u>	10 mg	189.00 USD
<u>D-3308</u>	<u>dextran, tetramethylrhodamine, 3000 MW, anionic, lysine fixable</u>	10 mg	214.00 USD
<u>D-3312</u>	<u>dextran, tetramethylrhodamine and biotin, 10,000 MW, lysine fixable (mini-ruby)</u>	10 mg	258.00 USD
<u>D-3328</u>	<u>dextran, Texas Red®, 3000 MW, lysine fixable</u>	10 mg	214.00 USD
<u>D-3329</u>	<u>dextran, Texas Red®, 3000 MW, neutral</u>	10 mg	189.00 USD
<u>D-3330</u>	<u>dextran, amino, 3000 MW</u>	100 mg	95.50 USD
<u>D34679</u>	<u>dextran, Alexa Fluor® 555, 10,000 MW, anionic, fixable</u>	5 mg	274.00 USD
<u>D34680</u>	<u>dextran, Alexa Fluor® 680, 10,000 MW, anionic, fixable</u>	5 mg	274.00 USD
<u>D34681</u>	<u>dextran, Alexa Fluor® 680, 3,000 MW, anionic</u>	2 mg	160.00 USD
<u>D34682</u>	<u>dextran, Alexa Fluor® 488, 3,000 MW, anionic</u>	2 mg	160.00 USD
<u>D-7132</u>	<u>dextran, Cascade Blue®, 3000 MW, anionic, lysine fixable</u>	10 mg	214.00 USD
<u>D-7135</u>	<u>dextran, biotin, 3000 MW, lysine fixable (BDA-3000)</u>	10 mg	214.00 USD
<u>D-7136</u>	<u>dextran, fluorescein, 500,000 MW, anionic, lysine fixable</u>	10 mg	189.00 USD
<u>D-7137</u>	<u>dextran, fluorescein, 2,000,000 MW, anionic, lysine fixable</u>	10 mg	189.00 USD
<u>D-7139</u>	<u>dextran, tetramethylrhodamine, 2,000,000 MW, lysine fixable</u>	10 mg	189.00 USD
<u>D-7142</u>	<u>dextran, biotin, 500,000 MW, lysine fixable (BDA-500,000)</u>	10 mg	189.00 USD
<u>D-7144</u>	<u>dextran, amino, 500,000 MW</u>	100 mg	95.50 USD
<u>D-7153</u>	<u>dextran, Rhodamine Green™, 10,000 MW, lysine fixable</u>	10 mg	258.00 USD
<u>D-7156</u>	<u>dextran, fluorescein and biotin, 3000 MW, anionic, lysine fixable (micro-emerald)</u>	5 mg	258.00 USD
<u>D-7162</u>	<u>dextran, tetramethylrhodamine and biotin, 3000 MW, lysine fixable (micro-ruby)</u>	5 mg	258.00 USD
<u>D-7163</u>	<u>dextran, Rhodamine Green™, 3000 MW</u>	5 mg	325.00 USD
<u>D-7168</u>	<u>dextran, BODIPY® FL, 10,000 MW, fixable</u>	5 mg	189.00 USD
<u>D-7170</u>	<u>dextran, Oregon Green® 488, 10,000 MW, anionic</u>	5 mg	189.00 USD
<u>D-7171</u>	<u>dextran, Oregon Green® 488, 10,000 MW, anionic, lysine fixable</u>	5 mg	214.00 USD
<u>D-7172</u>	<u>dextran, Oregon Green® 488, 70,000 MW, anionic</u>	5 mg	189.00 USD
<u>D-7173</u>	<u>dextran, Oregon Green® 488, 70,000 MW, anionic, lysine fixable</u>	5 mg	214.00 USD
<u>D-7176</u>	<u>dextran, Oregon Green® 514, 70,000 MW, anionic</u>	5 mg	189.00 USD
<u>D-7178</u>	<u>dextran, fluorescein and biotin, 10,000 MW, anionic, lysine fixable (mini-emerald)</u>	10 mg	257.00 USD
<u>N-7167</u>	<u>NeuroTrace® BDA-10,000 Neuronal Tracer Kit</u>	1 kit	259.00 USD
<u>P10361</u>	<u>dextran, pHrodo™ 10,000 MW for endocytosis</u>	0.5 mg	265.00 USD
<u>V-22916</u>	<u>Vybrant® Cell Lineage Tracing Kit</u>	1 kit	141.00 USD